




Liters/Min

Personal Best®
FULL RANGE
PEAK FLOW METER

**WITH INTEGRATED
ASTHMA MANAGEMENT
ZONE SYSTEM**

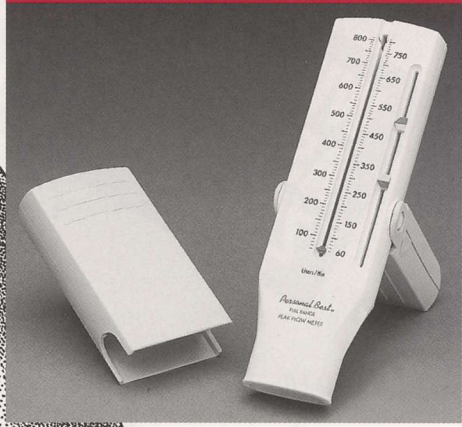

**IT'S NOT JUST
OUR NAME.
IT'S OUR PROMISE.**

PULMODATA
d.o.o.

Profesionalna medicinska oprema
razvoj, inženiring, zastopanje

Balantičeva 9
1234 MENGEŠ, SLOVENIJA
tel.: 01/723 95 11, fax: 01/723 95 15
e-mail: info@pulmodata.si
www.pulmodata.si

Personal Best®

Peak Flow Meter
with Integrated Asthma Management
Zone System

Product Specifications

Personal Best® meets or exceeds National Asthma Education and Prevention Program Technical Standards for Peak Flow Meters, based on American Thoracic Society Standardization of Spirometry 1994 Update (26 Waveforms).

Performance Data

Personal Best Standards	(whichever is greater)
Accuracy	± 10% or 20 L/min
Reproducibility	≤ 5% or 10 L/min
Interdevice Variability	≤ 10% or 20 L/min

Dimensions

Length	6.5"
Width	2.0"
Height	0.8"
Weight	3.0 oz

Materials

Main Body	Impact-resistant ABS plastic
Cover and Handle	High-density polypropylene
Scale	Hot-stamped, alcohol-resistant

Calibration Data

	Full Range	Low Range
Measurement Range	60–810 L/min	50–390 L/min
Resolution	10 L/min	5 L/min

Cleaning

Personal Best® can be washed with warm water and mild liquid soap or on the top rack of a dishwasher and air dried. Never boil or autoclave.

Life Span Data

Accuracy and reproducibility assured for a minimum of two years (±1%).



Respironics HealthScan
41 Canfield Road
Cedar Grove, New Jersey
USA • 07009-1201
1-800-962-1266
Telephone: 973-857-3414
Fax: 973-239-0831

Respironics International, Inc.
77, rue de Paris
92100 Boulogne Billancourt
FRANCE
Telephone: 33-01-55-60-19-80
Fax: 33-01-55-60-19-89

www.respironics.com



Quality Assurance

100% of Personal Best units are quality-control inspected at the factory before shipment.

Warranty

Personal Best is backed by a one-year unconditional warranty. Any unit performing unsatisfactorily for any reason during this period will be replaced free of charge by Respironics.

Three-Zone Asthma Management

The integrated three-zone asthma management system recommended by NAEPP helps patients comply with treatment regimens. The color-coded indicators (located to the right of the scale) can be adjusted to delineate a patient's green, yellow and red zones based on personal best peak flow. (These indicators have been designed to be difficult to move inadvertently.)

How Supplied

Each Personal Best unit is wrapped in a tamper-proof plastic bag with a patient record chart and instructions for use. Reorder Nos. 755 and 756 are packaged in a unit box and contain a Companion Card.

Ordering Information

Reorder No.	Product Description	Units/Case
HS755-012	Personal Best Full Range	12
HS756-012	Personal Best Low Range	12
HS755EU-010	Personal Best Full Range (CE Marked with labeling in 8 languages)	10
HS756EU-010	Personal Best Low Range (CE Marked with labeling in 8 languages)	10
HS714-200	Single-Use One-Way Valve Safety Mouthpieces for use with HealthScan Peak Flow Meters	200
HS757-100	Single-Use Mouthpieces for use with HealthScan Peak Flow Meters with Oval Mouthpieces	100

This product should be used under the supervision of a physician or licensed healthcare professional.

1003615 REV1

©2001 Respironics HealthScan, Inc.